

DECLARATION OF BLOOD PRESSURE MEASURING DEVICE EQUIVALENCE 2013

A SIGNED COPY WILL BE POSTED ON THE www.dableducational.org WEBSITE

SECTION A - Please complete all items.

I **Jacqueline Hart**,
Name of a Company Director

a Director of **Omron Healthcare Europe B.V.**,
Company name

hereby state that there are no differences that will affect blood pressure measuring accuracy between the

Maker^a **Omron Healthcare Co., Ltd.** Address **53, Kunotsubo, Terado-cho, Muko, Kyoto 617-0002 Japan**
 Manufacturer^b **Omron Healthcare Co., Ltd.** Address **53, Kunotsubo, Terado-cho, Muko, Kyoto 617-0002 Japan**
 Brand^c **Omron** Model^d **M6 Comfort IT (HEM-7322U-E)**

Blood pressure measuring device for which validation is claimed. If alternative model names are used, include all.

blood pressure measuring device and the validated blood pressure measuring device

Maker^a **Omron Healthcare Co., Ltd.** Address **53, Kunotsubo, Terado-cho, Muko, Kyoto 617-0002 Japan**
 Manufacturer^b **Omron Healthcare Co., Ltd.** Address **53, Kunotsubo, Terado-cho, Muko, Kyoto 617-0002 Japan**
 Brand^c **Omron** Model^d **M6 Comfort (HEM-7321-E)**

Existing validated blood pressure measuring device.

which has previously passed the **ESH2010** protocol, the results of which were published as follows:

dablEducational Trust; 2014 Jan 22. 4 p. Available from: <http://www.dableducational.org/Publications/2014/ESH-IP>
 2010 Validation of Omron M6 Comfort (HEM-7321-E).pdf

Full reference

The only differences between the devices involve the following components:

Tick one box for each item 1–18.

Part I	1	Algorithm for Oscillometric Measurements	Yes <input type="checkbox"/>	No <input checked="" type="checkbox"/>	N/A ^e <input type="checkbox"/>
	2	Algorithm for Auscultatory Measurements	Yes <input type="checkbox"/>	No <input type="checkbox"/>	N/A ^f <input checked="" type="checkbox"/>
	3	Artefact/Error Detection	Yes <input type="checkbox"/>	No <input checked="" type="checkbox"/>	
	4	Microphone(s)	Yes <input type="checkbox"/>	No <input type="checkbox"/>	N/A ^f <input checked="" type="checkbox"/>
	5	Pressure Transducer	Yes <input type="checkbox"/>	No <input checked="" type="checkbox"/>	
	6	Cuffs or Bladders	Yes <input type="checkbox"/>	No <input checked="" type="checkbox"/>	
	7	Inflation Mechanism	Yes <input type="checkbox"/>	No <input checked="" type="checkbox"/>	
	8	Deflation Mechanism	Yes <input type="checkbox"/>	No <input checked="" type="checkbox"/>	
Part II	9	Model Name or Number	Yes <input checked="" type="checkbox"/>	No <input type="checkbox"/>	
	10	Casing	Yes <input checked="" type="checkbox"/>	No <input type="checkbox"/>	
	11	Display	Yes <input checked="" type="checkbox"/>	No <input type="checkbox"/>	
	12	Carrying/Mounting Facilities	Yes <input type="checkbox"/>	No <input checked="" type="checkbox"/>	
	13	Software other than Algorithm	Yes <input checked="" type="checkbox"/>	No <input type="checkbox"/>	
	14	Memory Capacity/Number of stored measurements	Yes <input type="checkbox"/>	No <input checked="" type="checkbox"/>	
	15	Printing Facilities	Yes <input type="checkbox"/>	No <input type="checkbox"/>	N/A ^g <input checked="" type="checkbox"/>
	16	Communication Facilities	Yes <input checked="" type="checkbox"/>	No <input type="checkbox"/>	N/A ^g <input type="checkbox"/>
	17	Power Supply	Yes <input type="checkbox"/>	No <input checked="" type="checkbox"/>	
	18	Other Facilities	Yes <input checked="" type="checkbox"/>	No <input type="checkbox"/>	N/A ^g <input type="checkbox"/>

An explanation of each item ticked "Yes" must be included in **Section B** or on a separate sheet.

- Notes:
- a Provide the name and address of the actual maker of the device.
 - b Provide the name and address of the legal manufacturer of the device, even if it is the same as that of the maker.
 - c Provide the name of the brand under which it is sold, even if it is the same as that of the manufacturer or maker.
 - d Provide the model name. If alternative or internal model names are used, include all. Each device must be uniquely identifiable.
 - e Only tick N/A (Not Applicable) if neither device measures blood pressure using the oscillometric method.
 - f Only tick N/A (Not Applicable) if neither device measures blood pressure using the auscultatory method.
 - g Only tick N/A (Not Applicable) if neither device provides printing, communication or other facilities, as appropriate.

SECTION B An explanation for each item, 1 to 18, ticked "Yes" in Section A must be provided here or in an attached document. All differences between the devices must be described.

- 9) The model number is changed to M6 Comfort IT (HEM-7322U-E) from M6 Comfort (HEM-7321-E).
- 10) The USB port is added to the M6 Comfort IT.
- 11) The OK symbol, the transfer indicator and the DATA/FULL symbol are added to the M6 Comfort IT.
- 13) The M6 Comfort IT can manage to blood pressure value by software "BiLink".
- 16) The USB interface to connect with personal computer are added.
- 18) The M3 IT can manage to blood pressure value by software "BiLink".

SECTION C Please check that the following are included with the application

- A manual for the validated device
 - A manual for the device for which equivalence is being sought
 - An image of the validated device
 - An image of the device for which equivalence is being sought
 - An image of the screen layout of validated device*
 - An image of the screen layout of the device for which equivalence is being sought*
- * Screen layouts shown complete, and without obscuring labels or lines, in manuals need not be included separately.


SECTION D Complete all items, bar signatures and seal, online and print. Sign and seal it then send the original to our address below. Please email a signed copy of this form, together with the manuals and images for both devices, to info@dableducational.org.

Signature of Director 

Company Stamp/Seal

Name JACQUELINE HART

Date 27/11/14

Signature of Witness 




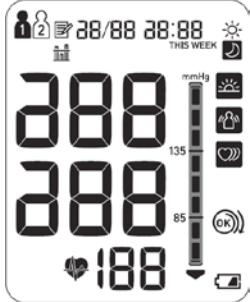


Name ANITA KEGSKES

Address Lusterburenstraat 146 1013GJ Amsterdam NL

OMRON HEALTHCARE EUROPE BV
 Scorpius 33
 NL-2132 LR Hoofddorp
 P.O.BOX 2050 NL-2130 GL Hoofddorp
 TEL +31-23 5544700
 FAX +31-23 5544701

OMRON HEALTHCARE EUROPE BV
 Scorpius 33
 NL-2132 LR Hoofddorp
 P.O.BOX 2050 NL-2130 GL Hoofddorp
 TEL +31-23 5544700
 FAX +31-23 5544701

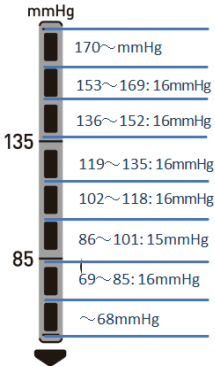
Comparison of the Omron M6 Comfort IT (HEM-7322U-E) with the Omron M6 Comfort (HEM-7321-E)

Devices	Omron M6 Comfort IT (HEM-7322U-E)	Omron M6 Comfort (HEM-7321-E)
Pictures		
Display		
Validation		ESH 2010
Device 1 Criteria	<p>Display/Symbols/Indicators</p> <p><i>Communication</i></p> <p>USB Connection ^{Query 3} 11, 16</p> <p>Data transmission completed  ^{Query 3} 11, 16</p> <p>Data awaiting transmission (Data/Full ) ^{Query 3} 11, 14, 16</p> <p>Algorithms</p> <p><i>Communication</i></p> <p>Data transfer to online database 13, 16</p> <p>Casing</p> <p><i>Ports</i></p> <p>USB port, cable and downloadable PC software 16, 18</p>	
Same Criteria	<p>Measurement</p> <p><i>Accuracy</i></p> <p>BP accuracy ± 3 mmHg 1, 5</p> <p>Pulse accuracy ± 5% 1, 5</p>	<p>Measurement</p> <p><i>Accuracy</i></p> <p>BP accuracy ± 3 mmHg 1, 5</p> <p>Pulse accuracy ± 5% 1, 5</p>

Devices	Omron M6 Comfort IT (HEM-7322U-E)	Omron M6 Comfort (HEM-7321-E)
Same Criteria (continued)	Measurement (continued)	Measurement (continued)
	<i>Method</i>	<i>Method</i>
	Oscillometric measurement method	Oscillometric measurement method
	BP 0 mmHg to 299 mmHg ^{Query 1}	BP 0 mmHg to 299 mmHg ^{Query 1}
	Pulse 40 bpm to 180 bpm	Pulse 40 bpm to 180 bpm
	Manually initiated measurements	Manually initiated measurements
	Measurements are from single inflations	Measurements are from single inflations
	Prevent storing of result (Guest mode)	Prevent storing of result (Guest mode)
	<i>Inflation</i>	<i>Inflation</i>
	Inflation 0 mmHg to 299 mmHg	Inflation 0 mmHg to 299 mmHg
	Automatic Inflation	Automatic Inflation
	Fuzzy Logic	Fuzzy Logic
	Press button if BP > 210 mmHg	Press button if BP > 210 mmHg
	<i>Deflation</i>	<i>Deflation</i>
	Automatic Deflation	Automatic Deflation
	<i>Cuffs</i>	<i>Cuffs</i>
	Large (Arm circ. 22 cm to 42 cm) No. HEM-FL31	Large (Arm circ. 22 cm to 42 cm) No. HEM-FL31
	<i>Measurement Records</i>	<i>Measurement Records</i>
	Memory: 100 measurements × 2 users (Guest not recorded)	Memory: 100 measurements × 2 users (Guest not recorded)
	Buttons/Switches	Buttons/Switches
	<i>Power</i>	<i>Power</i>
	On/Off with Start/Stop (Start/Stop Label)	On/Off with Start/Stop (Start/Stop Label)
	<i>Measurement Records</i>	<i>Measurement Records</i>
	Memory	Memory
	User ID slider	User ID slider
	<i>Function</i>	<i>Function</i>
	Date/Time set	Date/Time set
	Up and down	Up and down
	<i>Analysis</i>	<i>Analysis</i>
	Average	Average
Display/Symbols/Indicators	Display/Symbols/Indicators	
<i>Measurement Procedure</i>	<i>Measurement Procedure</i>	
Deflation symbol	Deflation symbol	
During Measurement: BP Level & Heartbeat	During Measurement: BP Level & Heartbeat	
<i>Post Measurement</i>	<i>Post Measurement</i>	
SBP, DBP and Pulse	SBP, DBP and Pulse	
Measurement error E1 E2 E3 E4 E5 Er	Measurement error E1 E2 E3 E4 E5 Er	

Devices	Omron M6 Comfort IT (HEM-7322U-E)	Omron M6 Comfort (HEM-7321-E)
Same Criteria (continued)	Display/Symbols/Indicators (continued)	Display/Symbols/Indicators (continued)
	<i>Post Measurement (continued)</i>	<i>Post Measurement (continued)</i>
	Hypertension (Indicator strip) ^{Query 2}	Hypertension (Indicator strip) ^{Query 2} 11, 13
	Hypertension (Indicator LEDs)	Hypertension (Indicator LEDs) 11, 13
	Morning hypertension	Morning hypertension 11, 13
	Memory zone average	Memory zone average 11, 13, 14
	7-day morning memory zone average	7-day morning memory zone average 11, 13, 14
	7-day evening memory zone average	7-day evening memory zone average 11, 13, 14
	Body movement error	Body movement error 3, 11, 13, 18
	Irregular heartbeat	Irregular heartbeat 11, 13, 18
	Correct cuff wrapping indicator (icon + 2 LEDs)	Correct cuff wrapping indicator (icon + 2 LEDs) 11, 13, 18
	<i>Measurement Records</i>	<i>Measurement Records</i>
	Memory icon	Memory icon 11
	Memory recall number (Replaces pulse rate momentarily)	Memory recall number (Replaces pulse rate momentarily) 11
	User (1, 2 and Guest)	User (1, 2 and Guest) 11
	<i>Date and Time</i>	<i>Date and Time</i>
	Date and Time	Date and Time 11
	Date and Time (During memory recall)	Date and Time (During memory recall) 11
	<i>Power</i>	<i>Power</i>
	Low & Exhausted battery	Low & Exhausted battery 11, 17
	Algorithms	Algorithms
	<i>Averages and Differences</i>	<i>Averages and Differences</i>
	Last 3 measurements (within 10 min) memory zone mean	Last 3 measurements (within 10 min) memory zone mean 13
	7-day morning and evening memory zone means (8 weeks)	7-day morning and evening memory zone means (8 weeks) 13
	<i>Diagnostic</i>	<i>Diagnostic</i>
	BP classification ^{Query 2}	BP classification ^{Query 2} 13
	135 / 85 mmHg thresholds	135 / 85 mmHg thresholds 13
	Morning hypertension	Morning hypertension 13
	Irregular heartbeat detection	Irregular heartbeat detection 13
	Body movement error detection	Body movement error detection 3, 13
<i>Functions</i>	<i>Functions</i>	
Correct cuff wrapping detection	Correct cuff wrapping detection 13	
Casing	Casing	
<i>Display</i>	<i>Display</i>	
Single screen display	Single screen display 10	
Segment LCD	Segment LCD 10	

Devices	Omron M6 Comfort IT (HEM-7322U-E)	Omron M6 Comfort (HEM-7321-E)
Same Criteria (continued)	<p>Casing (continued)</p> <p><i>Power</i></p> <p>4 “AA” batteries ~ 1000 measurements 17</p> <p>AC adapter (S-9515336-9 or UK-9983666-5) (Optional) 17</p> <p>Automatic switch-off when not used for 2 min 17</p> <p>Rechargeable batteries not permitted 17</p>	<p>Casing (continued)</p> <p><i>Power</i></p> <p>4 “AA” batteries ~ 1000 measurements 17</p> <p>AC adapter (S-9515336-9 or UK-9983666-5) (Optional) 17</p> <p>Automatic switch-off when not used for 2 min 17</p> <p>Rechargeable batteries not permitted 17</p>
Comparable Criteria		
Device 2 Criteria		

<p>Queries</p>	<p>1</p>	<p>Query Each of the manuals states that the blood pressure measurement range is 0 mmHg to 299 mmHg. They also state that the monitor will not inflate above 299 mmHg. This means that the measurement ranges must be below this.</p> <p>According to ISO 80601-2-30 (2009), the device must be capable (in non-neonatal mode) of indicating at least 60 mmHg to 230 mmHg for SBP and 40 mmHg to 130 mmHg for DBP (201.12.1.103), so specifying these separately is necessary. It also requires that the pressure ranges provided are the rated pressures (201.7.9.2.9 h and 201.12.1.101) and that those measurements that are outside of these ranges trigger a technical alarm (201.12.1.106).</p> <p>a) What are the SBP and DBP rated ranges for each of the devices?</p> <p>b) Are there technical alarm ranges for each of the devices?</p> <p>Response a) There is no SBP and DBP rated ranges because we have not defined the rated range of cuff pressure which is actually limited by measurement range of the pressure (not blood pressure) 0 to 299 mmHg. The capability to measure the required SBP and DBP range (201.12.1.103) are confirmed by technical validation test.</p> <p>b) There is no technical alarm because of the reason above.</p> <p>Query The rated ranges for SBP and DBP are simply the ranges, within the inflation range, for which SBP and DBP values are displayed. Where a pulse is detected close to the maximum inflation pressure or the wave envelope suggests SBP as being close to the maximum inflation pressure, it may be rejected an unreliable estimate of SBP. DBP estimates close to zero can be similarly rejected. These are the technical alarm conditions.</p> <p>The reply suggests that there are no upper or lower limits to either SBP or DBP i.e. 299 mmHg ≥ SBP > DBP ≥ 0 mmHg. Is this correct?</p> <p>Response Yes, this is correct.</p> <p>Comment The reply answers the query fully.</p>
	<p>2</p>	<p>Query Both devices have vertical indicator strips. Each strip contains 8 sections with a 135 mmHg mark between the 3rd and 4th sections and an 85 mmHg mark between the 6th and 7th sections. From the diagrams in the respective manuals, their function appears to be to display the pulse pressure, with a series of sections lit from one indicating a range of SBP values to one indicating a range of DBP values. However, no explanation is provided in either manual. What range of pressure does each section represent and are they the same for both devices?</p> <p>Response Please find the range of the pressure for each section as shown. These are same for both devices.</p> <p>Comment The reply answers the query fully.</p> 

	3	<p>Query The M6 Comfort IT (HEM-7322U-E) has Transfer, OK and Data/Full symbols that are not described in the manual. Please outline the functions of these symbols.</p> <p>Response The Transfer indicator \rightleftharpoons is displayed when the USB cable is connected.</p> <p>The OK symbol is displayed when the data transfer has finished.</p> <p>The DATA/Full symbol \blacksquare is not displayed, for the selected user, if there are fewer than 80 measurements not transferred. It will blink if there are between 80 and 99 measurements not transferred. It will be shown as a steady symbol if the memory is full (100 measurements) and none of them are transferred. <i>(M6 Comfort IT)</i></p> <p>Comment The reply answers the query fully.</p>
Note	1	The only difference between these devices is that the M6 Comfort IT (HEM-7322U-E) has additional features to allow measurements to be transferred to an online database.
Recommendation	Equivalence is Recommended	
Date	20 March 2014	